

Otto Bock®

QUALITY FOR LIFE




BRAVO*racer*

Find
out!



BRAVO *racer*

Find
out!

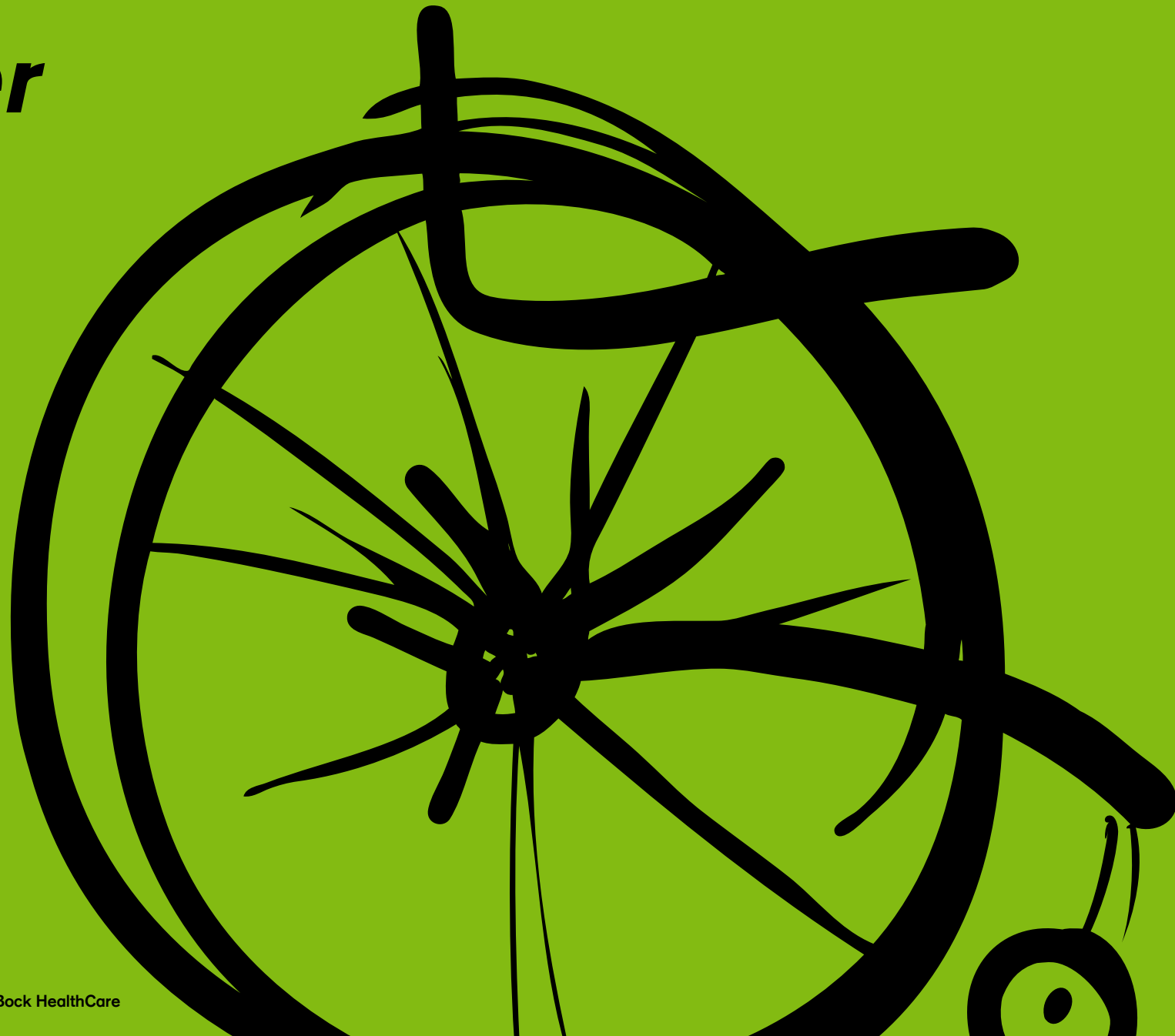



BRAVO!

The name promises a lot – and will keep its promises!

Who are we talking about? An active wheelchair for children and youths attracting attention everywhere with its fresh appearance. A companion with a stylish design that reflects a child's personality through its large selection of outfits. A practical, compact all-rounder that handles any everyday situation effortlessly and is easy to stow for transportation. An active wheelchair developed based on the body and growth data of children. It adapts to a child's requirements and implements the therapist's recommendations in a straightforward manner.

A dynamic-playful rigid-frame wheelchair now entering the large stage of children's wheelchairs – BRAVO*rac*er.

1
BRAVO!

BRAVO *racer*

Find
out!





2

BRAVOracer





2

BRAVOracer



Wheel lock

The non-slip cover assures easy, safe operation. Recessed in the side panel, the wheel lock lever is easy to reach but never in the way due to its compact size.



Anterior frame

The shapely, curved frame is more than just a fashion statement; it offers plenty of room for a broad footrest and effective leg abduction.



Back angle

The back rest is easy to unlock and fold down onto the seat, turning the wheelchair into a small, handy package.



Details




Clothing protector

The elegant side panel design provides additional support, protects clothing against dirt, and automatically retracts when the back rest is folded down. The concealed adjustment mechanism permits optimal adaptation to the user.



RightRun

The high-quality Schwalbe wheelchair tires run very easily and smoothly. These high-pressure air-filled tires are also available in a choice of colors!

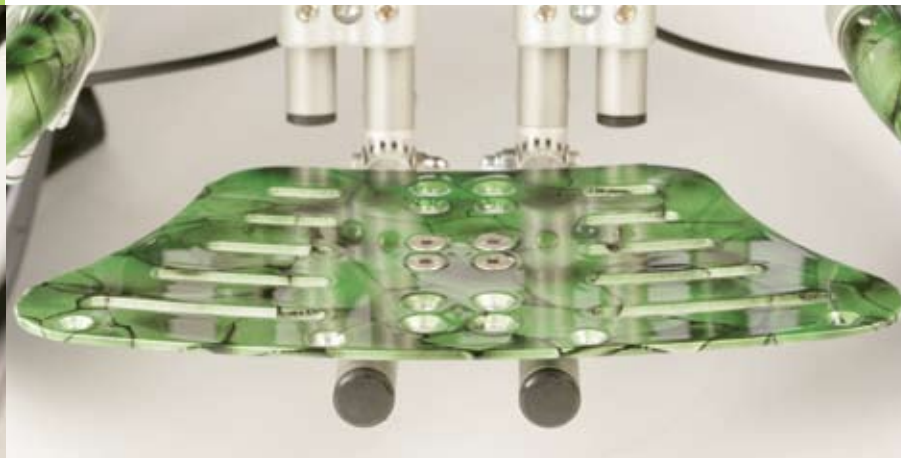

Footrest

Sturdy, generous, and easy to adjust for secure foot positioning. It can also be used as a step to help little ones get in or out!

2

BRAVO*racer*





2

BRAVO*racer*



Anti-tipper

The anti-tipper can be folded in or out in one simple step. It provides additional security for inexperienced drivers; in the passive position, it is hidden by the radius of the rear wheels.



Spoke protector

Useful and attractive: The spoke protector effectively prevents fingers from getting caught in the spokes and is an eye-catching accessory. It is shown here in the special „Gecko“ design.



Push bar

As an alternative to our selection of push handles, the push bar makes it easier to handle smaller wheelchairs in particular. When it is not needed, it can be removed in a few simple steps.



Options



Luminous casters

These smooth-running steering casters attract a lot of attention with their self-charging LEDs. As a visual gimmick or as an additional passive safety feature, the luminous casters are always a good choice.



Design clothing protector


This version of the side panel made from see-through materials offers additional color choices; it is supplied to match the wheelchair frame color.

2


BRAVO*racer*






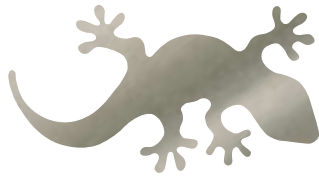


Brilliant Silver

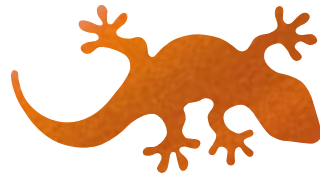


Sparkle Light Orange

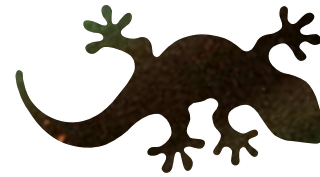


Snake Clear Red

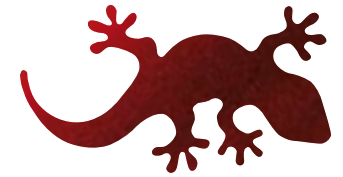


Candy Red

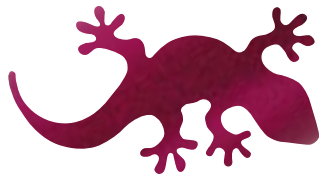

Dormant Vinho Sparkle

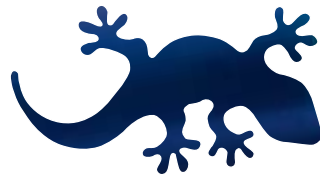

Dormant Marine Blue

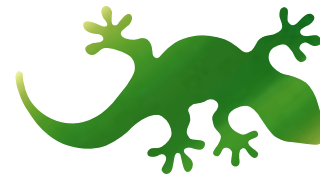

Dormant Apple Green



Clear Reptile Glossy



Colors



Spooky Yellow

Spooky Blue



3

Technical data





3

Technical data





Technical data

Seat width	20 cm – 36 cm
Seat depth	24 cm – 38 cm
Back height	20 cm – 40 cm
Back angle	-15° – +15° adjustable, can additionally be engaged in position +15°
Lower leg length	14,5 – 44,5
Seat height	31,5 – 46 (with 20" rear wheel); 34 – 48,5 (with 22" rear wheel); 36 – 51 (with 24" rear wheel)
Overall width	47 – 63; 48 – 66 (with 9° camber)
Overall length	61,5 – 76,5 (with 20" rear wheel); 64 – 79 (with 22" rear wheel); 66,5 – 81,5 (with 24" rear wheel)
Overall height (without push handles)	51,5 – 86 (with 20" rear wheel); 54 – 88,5 (with 22" rear wheel); 56,5 – 91 (with 24" rear wheel)
Weight	approx. 8,5 kg
Load capacity	60 kg

3

Technical data

BRAVOracer



Otto Bock®

QUALITY FOR LIFE

Otto Bock HealthCare GmbH

Max-Näder-Straße 15 · 37115 Duderstadt/Germany · Phone +49 5527 848-1590 · Fax +49 5527 848-1676
reha-export@ottobock.de · www.ottobock.com