

EasyStand® Magician® Comfy



The Way Kids Like to Stand



The Magician Comfy is a sit to stand stander for kids from 3'- 4'6" (91-137 cm) and up to 100 lbs (45 kg). It features Comfy-style wrap-around upholstery that provides the child with advanced support while standing and sitting.

To use the Magician Comfy, transfer the child from their wheelchair into the seat of the stander. After securing the positioning straps, kneepads and table, push the foot button to lift the child to the desired standing position using the gas cylinder lift.

The sit to stand design provides a supportive transition to standing, reducing anxiety for the child. Even when a child is unable to stand fully upright, they can safely stand partially upright, while being fully supported. The Magician Comfy can be adjusted to fit many different kids in multi-user environments, or can grow with your child at home.



Shown with options

John Cayden Murphy
Developmental Delay
47" (119 cm)

Magician Comfy Specifications

Height Range	3'-4'6"	91-137cm
Weight Limit	100 lbs.	45kg
Seat Depth Range		
-from seat pivot	10.5"-15.5"	27-39cm
Seat/Transfer Height		
-from floor	19"	48cm
Seat to Foot Plate Range	9.5"-16"	24-41cm
Table/Chest Pad Depth Range		
-from seat pivot	5"aft-6"forward	13-15cm
Table Height Range		
-from seat pivot	11"-22"	28-56cm
Chest Pad Size	7.5"x6"	19x15cm
Kneepad Depth Range		
-from front of seat	4"-10"	9-25cm
Knee Width Range	4.5"-7.5"	11-19cm
Seat Width Range	6"-13"	15-33cm
Back Height from Seat	11"-15"	28-38cm
Frame Color	Black Powder Coat	
Standard Upholstery Color	Royal Blue	
Seat Knee & Back Upholstery	Nylon Neoprene®	
Standard Desk Table Size	20"x22.5"	51x57cm
Standard Footprint	22"x33"	56x84cm
Weight of Basic Unit	95 lbs.	43kg

"Great for the classroom because it can be used as a desk and stander while providing great positioning. Gets the students out of their wheelchairs, provides them with the opportunity to weight bear and gives them confidence in an upright position. It's highly adjustable & can grow with the child."

Tiffany B., Physical Therapist - Greenfield, IN, USA



Standard Features

#7000F EasyStand Magician Comfy


The Magician Comfy comes standard with royal blue or ruby red wrap-around seat and back support with built-in seat belt and chest strap, kneepads, foot plates, 4" locking casters, clear acrylic table with front pad and gas spring lift.



Tool-free Adjustments



Seat Depth

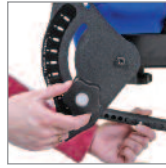

Back Angle



Back Height



Table/Front Pad



Kneepads



Foot Plates

Optional Accessories



#P7500M Mobile Option

Allows the child to self propel in the standing position. Includes locking push rims, belt drive and two front wheels. *Table size is 15" x 20" (38 x 51 cm); Mobile footprint is 23" x 35" (58 x 89 cm)*



#P80636 Foot Straps (Pair)

Hold the feet in place with Velcro® and d-ring adjustment. *Straps fit foot height of 2.5" - 4" (6-10 cm)*



#P80615 Secure Foot Option (Pair)

Velcro® and d-ring adjustment uses a four-point tightening system to hold the feet securely in place. *Straps fit foot height of 2.5" - 4" (6-10 cm)*



#P81071 Independent Kneepads (Pair)

Accommodate knee contractures and leg length discrepancies. *Kneepad size is 5" x 5" (13 x 13 cm); width range center to center is 6" - 12" (15 - 30 cm)*



#P5024 Table Swing-away Sleeve

Makes assisted transfers easier by moving the table away from the child.

Not available with the mobile option



#P80998 Angle Adjustable Table

Attaches to the existing table. Supports the child while they work at various angles from 10° to 45°.



#P81149 Chest Vest - Small

#P81150 Chest Vest - Medium

Neoprene® with camlock adjustment holds the child securely in position. *Small size is 9.5" L x 9" W (24 cm L x 23 cm W); medium size is 11" L x 9.5" W (28 cm L x 24 cm W)*



#P5070 Ratchet Handles

Can be used in place of existing knobs for ratchet motion and leverage. *Four small handles for table, two large handles for foot plates*



#P81315 Form-to-fit Upholstery (Royal Blue)

#P81316 Form-to-fit Upholstery (Ruby Red)

Seat and back include a formable metal insert for a customized fit.



#P80809 Tool Pouch

Attaches under the seat and holds the tools provided with the unit. *Tools not included*



#P81350 Head Support (Royal Blue)

#P81351 Head Support (Ruby Red)

Works as a prompt to limit head movement. *Height range from seat surface is 15" - 22" (33 - 56 cm); depth range from back surface 4" forward - 1" aft (20 cm forward - 3 cm aft)*



Headrest Ball Adapter

#P80869 Stealth® (.656" ball)

#P80870 Whitmyer (.75" ball)

#P80871 Whitmyer (.625" ball)

#P80872 Otto Bock™ (1" ball)